


SPORT-ELEC®

BODY BEAUTYFUL


BB2M

Mode d'emploi
User manual

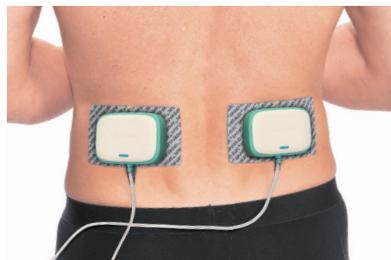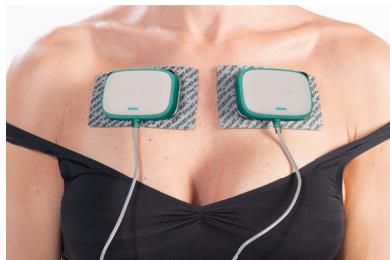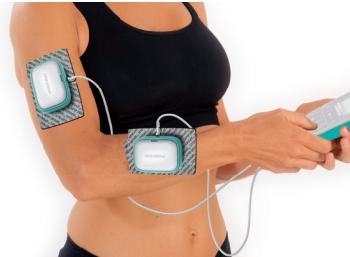DESCRIPTION DE L'APPAREIL / DEVICE DESCRIPTION

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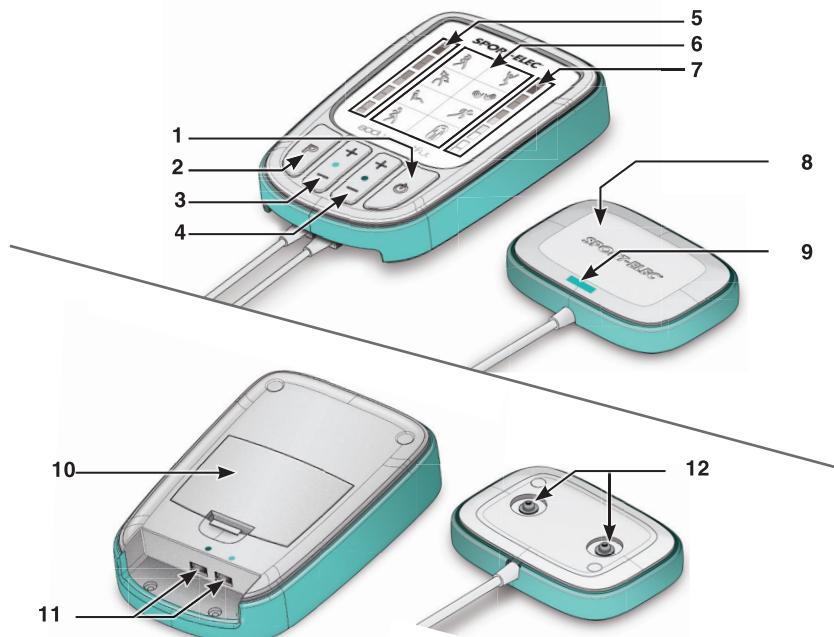
BODY BEAUTIFUL (BB2M)

DESCRIPTION OF THE UNIT (SEE PAGE 3)

1. «» ON/OFF key / Emergency stop.
2. «» Key Programs.
3. Buttons to adjust the power of A channel :
+ : to increase the power (from 0 to 72).
- : to decrease the power (from 72 to 0).
4. Buttons to adjust the power of B channel :
+ : to increase the power (from 0 to 72).
- : to decrease the power (from 72 to 0).
5. Indicator of A channel level power
6. Indicator of program
7. Indicator of B channel level power
8. Connexion module
9. Colour module (or channel) indicator
10. Compartment for the batteries
11. Connection sockets to plug the 2 modules (USB)
12. Fastener buttons

DESCRIPTION OF THE SCREEN

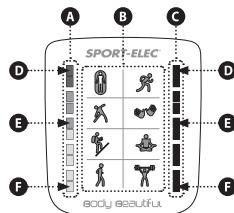- A** Indicator of A channel level power
- B** Indicator of program
- C** Indicator of B channel level power
- D** Used batteries : two red LED blink. Your device stops, even if a program is running. Change the batteries.
- E** Weak batteries indicator : two orange LED blink. Your device is still on operation.
- F** Disconnection warning light : blinks when the electrodes are disconnected.

BEFORE USING THIS UNIT

You have just acquired an effective and powerful professional unit. Take your time to discover how to use it.

- 1 Read the instructions carefully.
- 2 Before use, insert batteries in your device.
- 3 Test your unit on your biceps using the quick test procedure.

NON-THERAPEUTIC DEVICE

OPERATION QUICK TEST

TESTING YOUR UNIT ON YOUR BICEPS

- Connect 1 electrode to your device (see Chapter III. A).
- Test the device with the PG02 program.
- Increase power on this channel gradually. Press 1 for both the bearing from 0 to 1, then 8 consecutive presses to the following power levels.
- You would begin to feel contractions, the device is working.
- This discovery procedure shows your capacity in using the unit's level of performance.

PROGRAMS

PROG 1 PREPARATION

20 minutes prior to sport to prepare and warm up your muscles. The frequencies vary from 1 to 13Hz to provide a genuine massage effect.

Duration : 20 min

PROG 2 DEVELOPMENT

This program allows to improve the assets on the level of muscular form.

Duration : 40 min

PROG 3 FITNESS

Use this program to maintain the level of muscle endurance, resistance, and volume acquired.

Duration : 20 min

PROG 4 SUDOLORIS

This program relieves the muscular tensions effectively.

Duration : 20 min

PROG 5 SILHOUETTE SCULPT

40 minutes of body building to firm up, but also to increase muscular volume after a period of inactivity. This programme is ideal for keeping the body firm and for re-shaping your silhouette. Frequencies vary from 30 to 60Hz alternated with frequencies at 10Hz and 3Hz for active recovery to build up muscle resistance.

Duration : 40 min

PROG 6 RECOVERY

Thanks to fast frequencies changes this program allows the muscles to recover faster through relaxing effects.

Duration : 20 min

PROG 7 RELAXING MASSAGE

20 minutes to relieve muscular fatigue and allow you to recover after physical effort. The frequencies vary from 40 to 100Hz to provide the state of wellbeing of a relaxing massage.

Duration : 20 min

PROG 8 MUSCLE TONIC

20 minutes of more strenuous weight lifting to work on target areas. Frequencies vary from 70 to 100Hz alternated with frequencies at 3Hz for active recovery to build up tonicity and muscular strength.

Duration : 20 min

I. INSERTION OF BATTERIES (SEE PAGE 3)

Check that the device is switched OFF.

Open the batteries compartment, insert three 1,5 V AAA LR03 batteries into the compartment, comply with + and - indications then to close the case.

Do not mix used and good batteries. For an optimal use of your device, it is advised to use good quality alkaline batteries.

Do not use rechargeable batteries.

II. STARTING THE UNIT

A. PRESENTATION OF THE UNIT

Your electromuscular device is an high-tech appliance, it meets the current requirements in the field of medical application.

The interface «man machine» was simplified to the maximum; it is realized from lightened pictograms, indicators and function key. Among the eight pictograms the one enlightened in green displays you the program chosen by the «P» key.

Two scales indicator (1 per channel) displays you 72 levels of power.

You will find also 4 control keys (see diagram page 3):

- 1 «O» key red - press shortly

- 1 «P» key to choose the programs - press shortly

- 2 double keys «+/-» to adjust the power of stimulation (press shortly).

A long press (1 second) on the key «-» returns the chosen power rapidly to 0.

B. USING THE UNIT

STEP 1 ON / OFF «O»

A short press on the «O» key puts the device under or except power.

Every time when the power is on, the device carries out an auto-test on 11 essential safety parameters.

During the auto-test all the indicators are lit during ½ second.

If an internal defect in the device is detected, all the pictograms flash in green and the device stops. Then, return the device to After Sales Service.

The auto-test does not take into account the detection of the used batteries. When batteries are used, the auto-test is not activated. A long press (of more than 3 seconds) involves the release of the auto-test. In this case stop the device, please, by pressing the « ⌂ » key and then make a short press on the « ⌂ » key.

STEP 2 - PROGRAM SELECTION KEY « P »

To select a program, make one or more short presses on the « P » key, the programs scroll.

The programs are announced by pictograms.

As soon as a program is selected, corresponding pictogram lights.

Once the program is chosen, you can start it by increasing the power on one of power keys « + ».

During the program the « P » key is inactive.

STEP 3 - END OF THE PROGRAM - STOP DURING THE PROGRAM

At the end of the session the device stops automatically.

To stop the device during the program, put the 2 power keys back on zero or press the « ⌂ » key.

STEP 4 - ADJUSTMENT OF THE POWER

The power is divided into 10 levels : 5 green, 3 orange and 2 red levels.

To change the level press the « + » key. Press 1 for both the bearing from 0 to 1, then 8 consecutive presses to the following power levels.

You have 72 levels to regulate the power from 0 to 72.

STEP 5 - BATTERIES

If the batteries are good, there is no blinking.

If the batteries are weak, two orange indicators blink. Your device still functions.

If the batteries are used, two red indicators blink. Your device stops, even in the course of a program. Change the batteries.

STEP 6 - DISCONNECTION

Your device is equipped with a detector of electrodes disconnection.

If the contact between the human body and a single electrode is absent, the power is stopped on the concerned channel. You are warned of a disconnection by blinking of green indicator during 30 seconds.

If none of the electrodes is in contact with the body, the power cannot be emitted.

III. GENERAL INSTRUCTIONS

A. USING THE DEVICE WITH THE ELECTRODES

This function allows you to use the adhesive electrodes by connecting them directly to the modules.

You can use the device with 1 or 2 modules.

GETTING STARTED

- ▶ Connect each module to the unit by matching the color indicators (photo A page 4).
- ▶ Fix the electrodes to the module(s) using the snap fasteners (photo B page 4)

POSITIONING OF THE ELECTRODES

- ▶ Remove the plastic backing sheet from the electrodes (photos C page 4). Moisten them by spreading a drop of water onto the adhesive gel using your fingertips and wait 2 minutes.
- ▶ Position the electrode(s) onto your body.

NB : The electrodes are for personal use. For hygiene reasons, they must always be used by the same person.

STARTING THE PROGRAM

Start your device by following the « Starting the device » procedure (chapter II).

Contractions should be firm but must never cause any pain. If the electrodes cause you any discomfort : reduce the power levels, move the electrode(s) and then gradually increase the power. Also ensure that the electrodes are sufficiently moistened.

Now that the settings have been optimized, let your device get to work for the full duration of the program.

END OF PROGRAM

Once the program has finished, switch off the device by pressing the « Ø » button, then disconnect the modules. Gently remove the electrodes from your skin.

Moisten them slightly with a drop of water and replace them carefully onto their transparent plastic backing sheet. Keep the adhesive electrodes in a plastic bag in a cool place (for example, a separate compartment in your refrigerator).

B. USING THE DEVICE WITH THE BELT

The belt has 2 pairs of snap fasteners onto which 2 wireless electrodes specially designed by SPORT-ELEC® are clipped. To obtain new electrodes, contact your usual retailer or order them from our website www.sport-elec.com.

WARNING !


**⚠ Never use the belt without attaching the 2 electrodes provided.
Only use electrodes supplied by SPORT-ELEC®.**

PRESENTATION OF THE BELT

1 = Snap fasteners

2 = Wire storage flap
(with a lopp fixing device)

3 = Partie auto-agrippante

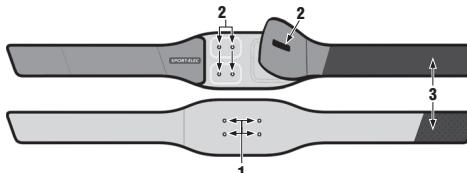CONNECTION AND POSITIONING OF THE BELT

- Check that your device is switched OFF.
- Open out the belt.
- Fix the electrodes to the inside of the belt using the snap fasteners.
- Place the belt around the waist and adjust it using the Velcro strip, so that there is a good contact between the electrodes and the skin. If required, use the belt extension.
- Once the belt is in position, fix 2 modules onto the snap fasteners.
- If you experience any tingling, switch off the device and check that the electrode is in contact with the skin. If the problem persists, contact the after-sales service : sav@sport-elec.com.

END OF PROGRAM

Once the program has finished, switch OFF the device and unplug it from the belt.

Remove the electrode and moisten them slightly with a drop of water before replacing them onto their transparent plastic backing sheets.

IV. PRECAUTIONS

- Do not use if pregnant, for 6 months after childbirth, in case of neuromuscular diseases or hemorrhagic, if a pacemaker.
- Hypersensitivity to electro-stimulation may cause a skin reaction in certain individuals.
- Keep out of the reach of children or mentally-impaired persons. This product contains small parts which can be easily swallowed.
- Only use BodyBeautiful with the accessories supplied and recommended by SPORT-ELEC® INSTITUT.
- Do not use BodyBeautiful when driving, operating machinery, or during any activity where electro-stimulation could result in the risk of injury to the user.
- Always position the belt in accordance with the information provided in this manual.
- Do not apply the electrodes to the lower part of your neck. Muscle contraction may be very strong and cause respiratory problems, as well as affect heart rate and blood pressure.
- Do not apply to the heart area, as electro-stimulation may cause heart rhythm disturbances.
- Do not place electrodes on the face, head or temples. The effects of electro-stimulation on the brain are unknown.
- Only apply the electrodes to clean and healthy skin. Do not carry out electro-stimulation on an open wound or on swollen, infected or inflamed areas, or on skin rashes (e.g. phlebitis, thrombophlebitis, varicose veins...).
- Do not use in the bath, shower or for any activity in contact with water.
- Do not use when sleeping.
- Do not use if you are fitted with an IUD contraceptive device, as contraceptive efficiency under electro-stimulation has not been established.
- Do not use on parts of the body other than those for which the accessories have been designed.

- The electrodes are for personal use. For hygiene reasons, they must always be used by the same person.
- Recycle or dispose of the battery in accordance with all federal, state and local laws. To avoid fire and explosion hazard, do not burn or incinerate the battery.
- The device is not suitable for use in the presence of flammable anesthetic mixtures with air or with oxygen or nitrous oxide.
- The device requires no calibration.
- The user must check that the equipment functions safely and see that it is in proper working condition before being used.
- Keep the equipment away from children and pets.

Special precautions should be taken in the following circumstances :

- Persons with suspected epilepsy.
- If you have a tendency to internal bleeding or bleeding following an injury or fracture.
- Following a surgical operation (electro-stimulation may disrupt the recovery procedure).
- In areas of the skin where sensitivity is impaired.

ADVERSE EFFECTS

- Skin irritation.
- Mild burns to the skin where the electrodes are applied.

MEDICAL ADVICE

Seek the advice of your doctor in the following cases :

- You suffer from serious back pain.
- You have a serious illness not mentioned in this manual.
- You have recently undergone a surgical procedure.
- You take insulin for diabetes.

RECOMMENDATIONS FOR USE

- You can restart a program cycle on another muscle zone.
- For optimum use, refer to the diagrams shown in this manual.
- The units can be cleaned using a cloth dampened with a disinfectant solution.

REPAIRS

If your device is not working correctly, stop using it and contact your retailer.

CONTACT SPORT-ELEC® INSTITUT OR AN AUTHORIZED DISTRIBUTOR IF :

- You experience any irritation, skin reaction, hypersensitivity or any other reaction. However, please note that any redness of the skin on contact with the electrodes is normal and disappears shortly after using your belt.
- Correct use of the electro-stimulation belt should not cause any discomfort. SPORT-ELEC® INSTITUT declines all responsibility if the user does not follow the instructions described in this manual.
- If for any reason you are in any doubt concerning the use of this device, consult your doctor prior to use.

RESTRICTIONS OF USE

- The belt must not be connected to another device or any other object other than the BodyBeautiful.
- The complete device must be used, with all detachable parts in place.
- Do not touch the inside of the belt when the device is switched on.
- Do not use the electro-stimulation belt when driving or operating machinery.
- Do not use in damp conditions.
- Use the device in a clean environment (dust and dirt-free...).
- Keep the device away from all sources of heat.
- To reposition the electro-stimulation belt during the session: stop the program in progress, loosen the belt, reposition it and restart once the belt has been readjusted.

WARNINGS

DO NOT USE THIS DEVICE

- If you carry a pacemaker.
- If you are pregnant.
- During first 6 weeks after baby birth.
- In case of neuromuscular diseases.
- In case of hemorrhagic diseases.

ASK FOR ADVICE FROM YOUR PHYSICIAN

- If you have electronic or metallic implants (pins...).
- In case of skin condition (wounds...).

MODIFYING THE UNIT IS FORBIDDEN

It is forbidden to modify the unit without the manufacturer's authorization.

If the unit is altered, a function check and an appropriate test must be carried out to ensure that it can still be used safely.

DO NOT APPLY ELECTRODES

- On the face, eyes, mouth or on the front of the neck.
- On one side to the other of (or through) the head.

POTENTIAL DANGERS

- The simultaneous connection of a patient to a highfrequency surgical unit may cause burns at thepoints of contact of the stimulator's electrodes, and the stimulator may be damaged.
- Operation in close proximity of a short wavetherapy unit may cause instability in the stimulator's output power.

V. TECHNICAL SPECIFICATIONS

UNIT

Class IIa medical device in accordance with Directive 93/42/CEE	
Type of current	Dual phase
Number of programmes	8 programs
Number of program variations (or phases)	55
Number of adjustable independent outputs	2
Current	7µA From 0 to 70 mA - From 0 to 1000 µA
Frequency range / Pulse width range	1 to 100 Hz / 160 to 320 µs
Power	3 x 1,5 V alkaline batteries AAA type LR03 (not provided)
Unit dimensions and weight	84 x 126 x 30 mm - 200 g
Automatic stop at end of each program	Yes - about 4 min 30 sec
Activation safety	100%
Absence of electrodes or of contact detection	Yes detection independent on each of the 2 outputs
Power indicator	Bar luminous graph
Weak and/or defective batteries indicator	Yes - luminous indicator
Absence contact electrodes Indicator	Yes - luminous indicator
Indicator of the programs	Luminous pictograms
Temperature of storage / Relative humidity	0 - 40 °C / 10 - 90 %
Operating temperature / Relative humidity	5 - 40 °C / 20 - 65 %
Connection device / electrodes	By wire (about 80 cm) with mini-B USB plug 2 snap fasteners per module (35 mm ←→ →)  Read the manuel thoroughly
Caution	2 years (unit only)
Guarantee	

BELT

References	CT-BB2M-EA
Size (unisex)	From 34 to 60
Length of the unfolded belt	122 cm
Fabric	35% sponge, 20% elastic fabric, 15% polyester, 15% nylon, 10% PVC + copper, 2% rubber, 3% nonwoven fabric

ADHESIVE ELECTRODES

Adhesive, hydrophilic and hypoallergenic electrodes	
Reference	EASF110x71
Contact area	67.5 x 47 mm (x 2)
Estimated life expectancy	40 uses
Please order new compatible adhesive electrodes	
on www.sport-elec.com .	

CE

EC declaration of conformity

Device in contact with the patient, excluding the heart area.

RoHS

These units are free of materials likely to have a negative impact on the environment.

IP 22

Protected against the ingress of foreign
≥ 12.5 mm and against dripping water when tilted up to 15°.

Caution : Read the instructions for use.

Maximum intensity per channel up to 1 kΩ : 90 mA.

Relative humidity for storage:
from 10 to 90%.

Temperature for storage:
from 0 to 40 °C.

Name of the manufacturer bringing the device to the market.

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The manufacturer respects the environmental law and contributes to the cost of recovery and recycling.

Follow the user manual instructions..

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This device is intended for domestic use.

Double insulated electrical appliance (>50V).

VI. GUARANTEE

We provide a 24 months guarantee for products supplied by us, from the date of purchase.

The adhesive electrodes are for personal use. For hygiene reasons, they can therefore neither be returned nor exchanged.

The life expectancy of the adhesive electrodes is estimated at approximately 40 uses.

Never attempt to repair or modify the device yourself or by a repairman not authorized by SPORT-ELEC®. The guarantee becomes null and void if repairs have been carried out by non-authorized persons or if the replacement parts used are not from the same manufacturer.

Similarly, only use compatible accessories manufactured by SPORT-ELEC®.

During the guarantee period we repair, free of charge, manufacturing defects or faulty parts, reserving the right to decide if certain parts need repairing or replacing, or if the device itself must be replaced.

Any damage resulting from improper use (knocks, plugging into a direct current, wrong voltage...) does not fall under the guarantee; normal wear and tear does not affect the proper functioning or value of the device.

The guarantee is only valid if the date of purchase and the stamp and signature of the retailer appear on the guarantee certificate or upon presentation of the till receipt (see « guarantee certificate » at the end of this manual).

No other claims fall under the guarantee, unless otherwise provided by law.

ACCESSORIES GUARANTEE

The electrodes are designed for personal use. For reasons of hygiene, there will be no trade-in nor exchange.

DO NOT FORGET TO ACTIVATE YOUR WARRANTY ONLINE: www.sport-elec.com/activation

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► GUIDANCE AND MANUFACTURER'S DECLARATION – ELECTROMAGNETIC EMISSIONS

The Body Beautiful is intended for use in the electromagnetic environment specified below. The customer or the user of the Body Beautiful should assure that it is used in such an environment.

EMISSIONS	COMPLIANCE	ELECTROMAGNETIC ENVIRONMENT - GUIDANCE
RF emissions CISPR 11	Group 1	The Body Beautiful uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class B	
Harmonic emissions IEC 61000-3-2	Conform Class A	
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Conform	The Body Beautiful is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.

► GUIDANCE AND MANUFACTURER'S DECLARATION
ELECTROMAGNETIC IMMUNITY

The Body Beautiful is intended for use in the electromagnetic environment specified below. The customer or the user of the Body Beautiful should assure that it is used in such an environment.

IMMUNITY TEST	IEC 60601 TEST LEVEL	COMPLIANCE LEVEL	ELECTROMAGNETIC ENVIRONMENT - GUIDANCE
Electrostatic discharge (ESD) IEC 61000-4-2	±6 kV contact ±8 kV air	±6 kV contact ±8 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30 %.
Electrical fast transient/burst IEC 61000-4-4	±2 kV for power supply lines	Not applicable	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	±1 kV line(s) and neutral	Not applicable	
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	<5 % Ur ($>95\%$ dip in Ur) for 0.5 cycle 40 % Ur (60 % dip in Ur) for 5 cycles 70 % Ur (30 % dip in Ur) for 25 cycles <5 % Ur ($>95\%$ dip in Ur) for 5s	Not applicable Not applicable Not applicable Not applicable	Mains power quality should be that of a typical commercial or hospital environment. If a dips or an interruption of mains power occurs, the current of the Body Beautiful may be dropped off from normal level, it may be necessary to use uninterruptible power supply or a battery
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	3 A/m	3A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.

NOTE - Ur is the a.c. mains voltage prior to application of the test level

► GUIDANCE AND MANUFACTURER'S DECLARATION
ELECTROMAGNETIC IMMUNITY

The Body Beautiful is intended for use in the electromagnetic environment specified below. The customer or the user of the Body Beautiful should assure that it is used in such an environment.

IMMUNITY TEST	IEC 60601 TEST LEVEL	COMPLIANCE LEVEL	ELECTROMAGNETIC ENVIRONMENT – GUIDANCE
Conducted RF IEC 61000-4-6	3 Vrms 150 kHz to 80 MHz	3 Vrms	<p>Portable and mobile RF communications equipment should be used no closer to any part of the Body Beautiful, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.</p> <p>Recommended separation distance $d=1.67\sqrt{P}$ $d=1.67/\sqrt{P}$ 80 MHz to 800 MHz $d=2.333/\sqrt{P}$ 800MHz to 2.5 GHz</p> <p>where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation Distance in metres (m).</p> <p>Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey, should be less than the compliance level in each frequency range.b Interference may occur in the vicinity of equipment marked with the following symbol:</p> 
Radiated RF IEC 61000-4-3	3 V/m 80 MHz to 2.5 GHz	3 V/m	

NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies.

NOTE 2 These guidelines may not apply in all situations.

Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

a. Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the Body Beautiful is used exceeds the applicable RF compliance level above, the Body Beautiful should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the Body Beautiful.

b. Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m.

► RECOMMENDED SEPARATION DISTANCES BETWEEN PORTABLE AND MOBILE RF COMMUNICATIONS EQUIPMENT AND THE FREE ACTION

The Body Beautiful is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of Body Beautiful can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the Body Beautiful as recommended below, according to the maximum output power of the communications equipment.

RATED MAXIMUM OUTPUT POWER OF TRANSMITTER (W)	SEPARATION DISTANCE ACCORDING TO FREQUENCY OF TRANSMITTER (M)		
	150 KHZ TO 80 MHZ $d = 1.17\sqrt{P}$	80 MHZ TO 800 MHZ $d = 1.17\sqrt{P}$	800 MHZ TO 2.5 GHZ $d = 2.33\sqrt{P}$
0.01	0.117	0.117	0.233
0.1	0.369	0.369	0.738
1	1.167	1.167	2.333
10	3.690	3.690	7.377
100	11.67	11.67	23.33

For transmitters rated at a maximum output power not listed above, the recommended separation distance d in metres (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

NOTE 1 At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people

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